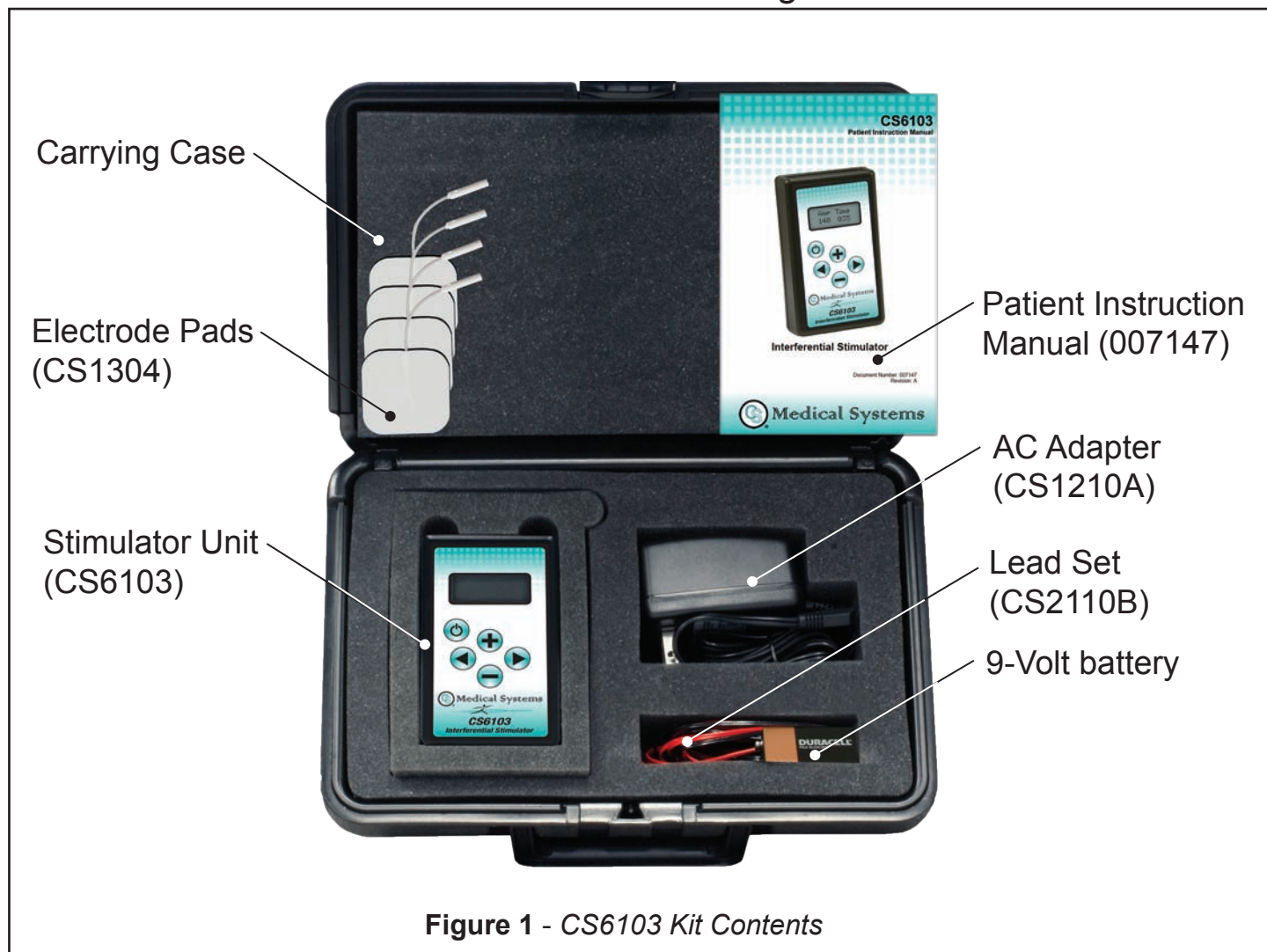


# CS6103 Interferential (IF) Stimulator Kit Quick Start Guide

The CS6103 kit includes the items shown in Figure 1.


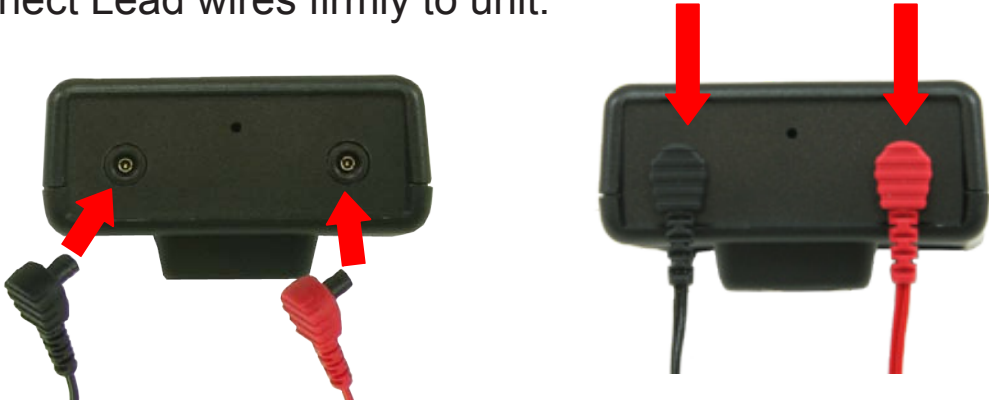
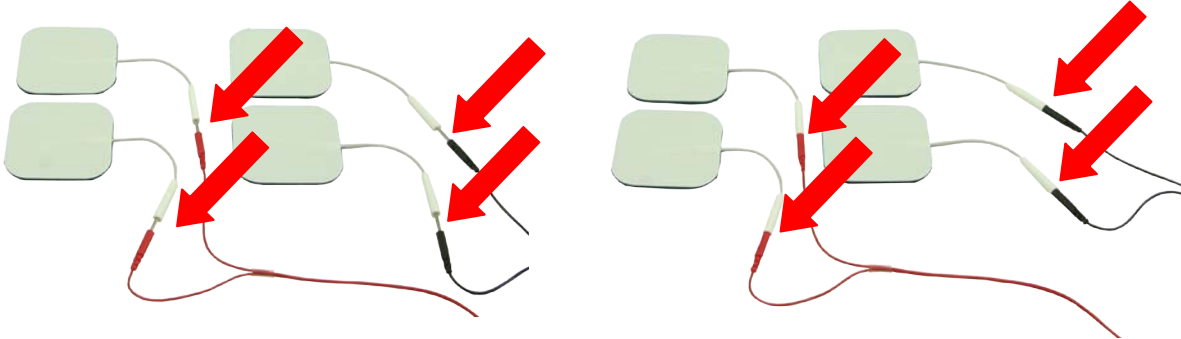








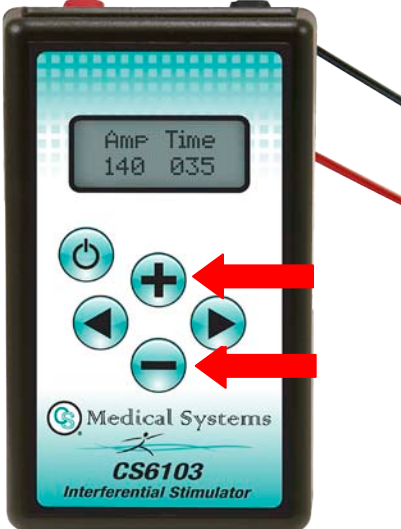
Once you have verified the kit contents, you can follow the steps to begin using your Interferential Stimulator.









Operating details can be found in the CS6103 Patient Instruction Manual, document number 007147, which is included as part of the kit.

Document Number: 007154  
Revision: A

Step	Action
<p>1</p>	<p>Provide power by:</p> <ul style="list-style-type: none"> <li>- Inserting a 9-volt battery, or</li> <li>- Plugging the AC Adapter into the unit and an AC receptacle.</li> </ul> 
<p>2</p>	<p>Connect Lead wires firmly to unit.</p> 
<p>3</p>	<p>Connect the Lead wires to Electrodes and place Electrodes per directions of healthcare provider.</p> 

Step	Action														
4	Press  to turn on stimulator. 														
5	Press  to move to Mode (Mod) field. 														
6	Press  or  to select the mode of operation prescribed by healthcare provider. Modes are: Std-36R-36A-Cnt-Stm-Mem. <table border="1" data-bbox="293 1507 1040 1822"> <thead> <tr> <th>Mode</th> <th>Treatment Time</th> </tr> </thead> <tbody> <tr> <td>Std - Standard</td> <td>20 minutes</td> </tr> <tr> <td>36R</td> <td>30 minutes</td> </tr> <tr> <td>36A</td> <td>60 minutes</td> </tr> <tr> <td>Cnt - Continuous</td> <td>60 minutes</td> </tr> <tr> <td>Stm - Muscle Stimulation</td> <td>20 minutes</td> </tr> <tr> <td>Mem - Memory</td> <td>Programmable</td> </tr> </tbody> </table> 	Mode	Treatment Time	Std - Standard	20 minutes	36R	30 minutes	36A	60 minutes	Cnt - Continuous	60 minutes	Stm - Muscle Stimulation	20 minutes	Mem - Memory	Programmable
Mode	Treatment Time														
Std - Standard	20 minutes														
36R	30 minutes														
36A	60 minutes														
Cnt - Continuous	60 minutes														
Stm - Muscle Stimulation	20 minutes														
Mem - Memory	Programmable														

Step	Action
7	Press  to move to the amplitude (Amp) field. 
8	Press  or  to set the Amp value prescribed by healthcare provider.  <div style="border: 1px solid black; padding: 5px; display: inline-block;">  If the Amp value is zero (0), the unit will not begin treatment.                     </div> 
9	Treatment will begin after the stimulator is idle for approximately 10 seconds.
10	The face of the stimulator will display the Amp value and count-down timer. <div style="border: 1px solid black; padding: 10px; margin: 10px auto; width: fit-content;"> <p style="text-align: center;">AMP    Time</p> <p style="text-align: center;">140    035</p> </div>
11	Treatment ends when timer reaches zero (0).
12	The unit will power off after 5 minutes.
13	When complete, carefully remove lead wires and electrodes.
14	Place the unit and all accessories in the provided carrying case.